

# Direct Heating Incubators




ZXDP-B2050S

## Descriptions:

ZXDP-series direct heating incubators offer an economical incubation method to a variety of micro-organism cultures. As a standard feature, all incubators have a high-tech microprocessor controller that maintains accurate temperature control (PID) of the chamber in the range of ambient +5°C to 65°C. With an accuracy of  $\pm 0.1^\circ\text{C}$ , a user-friendly sealed control panel allows easy digital setting and fast readout of time and temperature as well as an alarm signal if there is a deviation from set point parameters. The processor of our incubators can store up to 9 programs-18 steps-into a non-volatile memory. Programs can be used for running multiple "ramp and soak" cycles up to 99 times. Of course the incubators can also be used as a constant (one) temperature incubator.

## Features:

- P.I.D. microprocessor ensures the precision of temperature control under both fixed value mode and program mode.
- Fully programmable through the keyboard (9 programs/18 steps).
- Forced-air circulation provides uniformity of  $\pm 0.5^\circ\text{C}@37^\circ\text{C}$ .
- Audible and visible alarm on over temperature.
- Password protection against unauthorized access to all settings.
- Non-volatile memory retains preset parameters after an accidental power interruption.
- Large LCD display for easy parameters setting and fast readout.
- Inner glass door provides complete visibility to the chamber.
- Standard 2 grids included
- Low noise level: 38dB Rating (A)

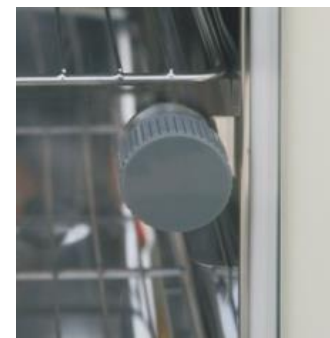

Inner Glass Door

### Specifications:

<b>Model</b>	<b>ZXDP-B2050S</b>
Heat Mode	Direct Heat (With Gentle Circulation)
Volume (L)	50
Temperature	ambient+ 5 to 65°C
Temperature Accuracy	± 0.1°C.
Temperature Uniformity	±0.5°C@37°C
Alarm	Enabled
Timer	0-999 minutes
Setting	Digital
Display	LCD
Grid Included	2 (max 4)
Grid Size(mm) (WxD)	330x345
Distance Between Grids (mm)	80
Inner Dimensions (mm) (WxDxH)	350x350x410
Exterior Dimensions (mm) (WxDxH)	470x520x785
Packing Dimensions (mm) (WxDxH)	540x590x945
Net Weight (kg)	33/57
Power (W)	200
Electricity	220-240V 50/60 Hz
Approval	CE, ISO

### Ordering Information:

Part NO.	Description	Part NO.	Description
ZXDP-B2050S	ZXDP-B2050S,50L, Direct Heat Incubator, A+5-65°C	P9001	Grid Plate for ZXDP-B2050, S/S
※ S/S: Stainless Steel			

Labwit quality assurance program demands the continuous assessment and improvement of all Labwit products. Information in this leaflet could thus change without notification and does not constitute a product specification. Please contact Labwit or their representatives if you require more details