

Direct Heating Incubators




ZXDP-B2080

Descriptions:

ZXDP-series direct heating incubators offer an economical incubation method to a variety of micro-organism cultures. As a standard feature, all incubators have a high-tech microprocessor controller that maintains accurate temperature control (PID) of the chamber in the range of ambient +5°C to 65°C. With an accuracy of $\pm 0.1^{\circ}\text{C}$, a user-friendly sealed control panel allows easy digital setting and fast readout of time and temperature as well as an alarm signal if there is a deviation from set point parameters. The processor of our incubators can store up to 9 programs-18 steps-into a non-volatile memory. Programs can be used for running multiple "ramp and soak" cycles up to 99 times. Of course the incubators can also be used as a constant (one) temperature incubator.

Features:

- P.I.D. microprocessor ensures the precision of temperature control under both fixed value mode and program mode.
- Fully programmable through the keyboard (9 programs/18 steps).
- Forced-air circulation provides uniformity of $\pm 0.5^{\circ}\text{C}@37^{\circ}\text{C}$.
- Audible and visible alarm on over temperature.
- Password protection against unauthorized access to all settings.
- Non-volatile memory retains preset parameters after an accidental power interruption.
- Large LCD display for easy parameters setting and fast readout.
- Inner glass door provides complete visibility to the chamber.
- Standard 2 grids included

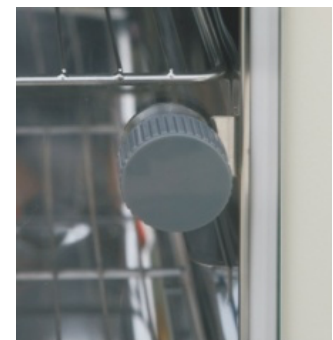

Inner Glass Door

Specifications:

Model	ZXDP-B2050	ZXDP-B2080	ZXDP-B2160	ZXDP-B2270
Heat Mode	Direct Heat (With Gentle Circulation)			
Volume (L)	50	80	160	270
Temperature	ambient+ 5 to 65°C			
Temperature Display Accuracy	± 0.1°C.	± 0.1°C.	± 0.1°C.	± 0.1°C.
Temperature Measuring Accuracy	± 1.0°C.	± 1.0°C.	± 1.0°C.	± 1.0°C.
Temperature Uniformity	±0.5°C@37°C	±0.5°C@37°C	±0.5°C@37°C	±0.5°C@37°C
Alarm	Enabled	Enabled	Enabled	Enabled
Timer	0-999 minutes	0-999 minutes	0-999 minutes	0-999 minutes
Setting	Digital	Digital	Digital	Digital
Display	LCD	LCD	LCD	LCD
Grid Included	2 (max 4)	2 (max 4)	2 (max 4)	2 (max 4)
Grid Size(mm) (WxD)	330x345	380x395	480x495	530x595
Distance Between Grids (mm)	80	100	130	160
Inner Dimensions (mm) (WxDxH)	350x350x410	400x400x500	500x500x650	600x550x820
Exterior Dimensions (mm) (WxDxH)	470x520x785	520x570x880	620x670x1030	740x740x1280
Packing Dimensions (mm) (WxDxH)	540x590x945	590x640x1040	690x740x1190	810x810x1440
Net Weight (kg)	33/57	40/73	57/94	90/130
Power (W)	200	250	380	550
Electricity	220-240V 50/60 Hz			
Approval	CE, ISO			

Ordering Information:

Part NO.	Description	Part NO.	Description
ZXDP-B2050	ZXDP-B2050,50L, Direct Heat Incubator, A+5-65°C	P9001	Grid Plate for ZXDP-B2050, S/S
ZXDP-B2080	ZXDP-B2080,80L, Direct Heat Incubator, A+5-65°C	P9002	Grid Plate for ZXDP-B2080, S/S
ZXDP-B2160	ZXDP-B2160,160L, Direct Heat Incubator, A+5-65°C	P9004	Grid Plate for ZXDP-B2160, S/S
ZXDP-B2270	ZXDP-B2270,270L, Direct Heat Incubator, A+5-65°C	P9005	Grid Plate for ZXDP-B2270, S/S

※ S/S: Stainless Steel

Labwit quality assurance program demands the continuous assessment and improvement of all Labwit products. Information in this leaflet could thus change without notification and does not constitute a product specification. Please contact Labwit or their representatives if you require more details